



Manual Large Volume Pipette Controller


Item Code: LIVPF100

The Manual Large Volume Pipette Controller may be used with glass or serological pipettes to facilitate the pipetting of a wide range of liquids.

Pipettes within the volume range of 0.1ml to 100ml can be used.


Special Features:

- *Ergonomically designed to allow better control of pipetting*
- *Made of durable, light weight material*
- *Wide volume range, 0.1 – 100ml*
- *Chemically resistant*

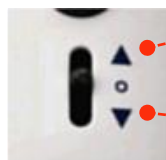

Filling – Squeeze the suction bellows and push up on the lever.

Dispensing – to dispense the fluid press down on the lever.

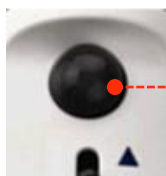

Blow-Out Function – For rapid expulsion of liquid, press the blow button.



Components

